

## Lactate measuring - what for?

- Whether for healthiness, fitness, recreation or match preparation: Only a targeted and healthy training leads to real success!
- The individual optimum training level under strain can precisely be determined only by measuring the increasing lactate value in blood:
- Inefficiency of training can be realized and avoided in time!
- Remarkable increase in endurance performance, optimum fat catabolism and a more balanced activity of the heart through training control!

## Products & accessories:

### Lactate Scout SOLO

Lactate Scout pocket device, 2x AAA/LR03 batteries, instruction manual

### Lactate Scout START

Lactate Scout pocket device, 2x AAA/LR03 batteries, carrying case with belt, lancing device, 25 lancets, 30 alcohol pads, 6 test strips, control solution 2,5ml and instruction manual with poster, service-CD

### Big Pack 24

24 test strips

### Optional:

- Big Pack 48+ test strips
- Big Pack 72+ test strips
- lancets
- lancing device
- alcohol pads
- control solution
- carrying case
- belt
- PC Packs

(All offers are subject to change and errors excepted, rights for alternative offers are reserved. All customers are bound by the terms and conditions of EKF-diagnostic sales GmbH)

**For further information please contact our sales team:**

Sales:

**EKF**  
DIAGNOSTIC

Sales:  
EKF - diagnostic sales GmbH  
Ebendorfer Chaussee 3  
D-39179 Barleben/Magdeburg  
Tel. +49(0)39203-785-0  
Fax +49(0)39203-785-24  
www.EKF-diagnostic.de


**SensLab**

producer:  
SensLab GmbH  
a Member of EKF - Group  
Bautzner Str. 67  
D - 04347 Leipzig  
Tel. +49(0)341-2341841  
Fax +49(0)341-2341840  
www.senslab.de



# LACTATE SCOUT

**Performance  
under control**

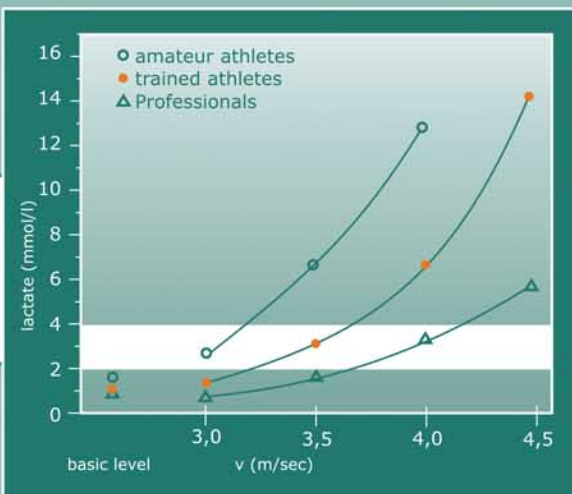




## Lactate SCOUT

- The companion for effective training!
- Fast & easy to operate with teststrips and comfortable scrollwheel
- Just a drop of blood (0.5 µl) is required
- Highly credible results within 15 sec.
- Measuring range: 0,5 - 25,0 mmol/l
- VC: 3 - 8% (dependent on concentration)
- Memory for 250 measurements, PC - link & individual training programme

## Fast & easy measuring:

- Just pick a test strip out of the vial & insert into the Lactate Scout
- Clean, sanitize & dry the earlobe or fingertip carefully
- Prick the area for blood withdrawal
- Hold the Lactate Scout with test strip onto the blood droplet for a moment
- Only 0,5 µl of blood will be absorbed, the measurement occurs automatically
- After 15 sec measuring duration, the result is displayed and stored
- Manifold options for gradual tests & Software analysis on the PC

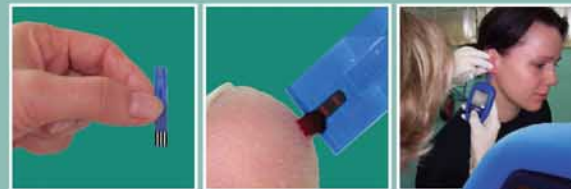

## Can be used anywhere:

- in sports medicine diagnosis for medical supervision
- in physiotherapy & recreation
- in sports institutions, fitness studios, for sports clubs & personal coaching
- in training of individuals & teams, recreational & competitive sports
- as a flexible & inexpensive complement to laboratory analysis
- can also be used privately, comfortable & hygienic in application
- accredited as a medical device according to EU Directive IVD (CE 0483)

# LACTATE SCOUT

- Professional Lactate measuring fast, accurate & easy!
- Metabolism based diagnostics for optimum endurance & fitness!
- The intelligent efficiency analysis for therapy and professional sport!
- A sports-medical companion for trainers, amateurs and clubs!

